



Euroseal® Valida

Everything you could desire from a sealing device. Plus some extras.

The new Euroseal® Valida is the thermosealing machine that “validates” the sealing cycle, memorising the sealing data and the result of operations to an SD card, as required by the EN11607-2 and EN868-5 Standards. Euroseal® Valida can be connected directly to the PC via Ethernet and WiFi system for instant transfer of sealing data. The innovative touch-screen display with dedicated menu allows the user to achieve simple and instant interfacing with the sealing machine and choose from the preset cycles and a free cycle, where it is possible to change the sealing temperature according to the type of roll and pouch being used.

- Sealing strip: 12 mm
- Sealing area width: 310 mm
- Weight: 8 kg
- Dimensions:
- standard L 500 x D 490 x H 315 mm
- with double roll holder L 500 x D 490 x H 420 mm
- with wall mounted roll holder L 500 x D 310 x H 210 mm



Euroseal® 2001 Plus

Electronic thermosealing machine for dentists

Electronic thermosealing machine at continuous cycle and constant temperature. The reduced overall dimensions allow an easy fitting in standard sterilization room cabinets. The metal roll-locking bar is designed to improve machine performance, holding the roll permanently still and inside the welding area. The smooth surface and ergonomic design allow easy machine cleaning.

- Sealing strip: 12 mm
- Sealing area width: 310 mm
- Weight: 7,6 kg
- Overall dimensions: L 445 x H 198 x D 385 mm



Euroseal® Check

Sterilization rolls for dental clinics

Bar heat thermosealing machine control test which allows to verify that rolls and pouches are correctly sealed, conforming to the EN ISO 11607-2.

- Box: 100 tests**
- Pack: carton of 40 boxes**



Euroseal®

Design, quality, simplicity and versatility in only one thermosealing machine.

Transparent polycarbonate cover to allow the possibility of carrying out the sealing procedures with better precision. Double roll holder for the simultaneous use of various sterilization rolls, hanging up to the wall. Adjustment of height and depth of the roll holder in order to adapt to the needs of each single dentist surgery. Integrated support panel to facilitate the sealing procedures and allows to work in a more functional manner. Anti roll return blocking system. Safety system to avoid the burning of the rolls. Visual and acoustic signals during the use. In full compliance with UNI EN 868-5 and EN 11607-2.

- Sealing strip: 12 mm
- Sealing area width: 310 mm
- Weight: 7,1 kg
- Dimensions:
- standard L 440 x D 390 x H 315 mm
- with double roll holder L 440 x D 390 x H 420 mm
- with wall mounted roll holder L 440 x D 300 x H 200 mm

STERILITY UP TO 6 MONTHS



Eurosteril® sterilization rolls and pouches

Design, quality, simplicity and versatility in only one thermosealing machine.

NEW paper rolls Medical high weight (60g/m²) and double layer of polyester film / polypropylene. The perfect solution to safely store your instruments. In full compliance with 93/42/EEC (and subsequent amendment 2007/47/EC) and EN 868-1.

- Measure available:**
- Rolls: 5,5-7,5-10-15-20-25-30-40 cm x 200 m and 20x5 cm x 100 m**
- Self-sealing pouches: 9x 25 cm; 14 x 26 cm; 19 x 33 cm**
- Pouches: 5,5x25 cm; 14 x 26 cm; 19 x 33 cm**



Eurosonic® 4D

Digital ultrasonic tank for dental clinics

Eurosonic 4D is Euronda Sweep Mode technology. Is a digital and electronic ultrasonic tank appaeared with three piezoceramic transducers. Its operating system is based on a highly user-friendly logic with luminous indicators on the controller and LED indicators linked to the various cycles, times and temperatures. Eurosonic® 4D is an electronic and digital ultrasonic tank with two operational modes, a temperature range of from 30°C to 60°C and time range of from 10 to 40 minutes.

- Large capacity stainless steel tank: 3,5 litres
- Overall dimensions: L 320 x H 240 x D 170 mm
- Weight: 4,5 kg
- Available accessories: perforated stainless steel basket, smoke-grey anti-deformation cover with non-drip system, 1 x 600 ml beaker holder, 3 x 600 ml beaker holder, glass beakers.



Eurosonic® 3D

Digital ultrasonic tank for dental clinics

Digital ultrasonic tank featuring a constant operating temperature of 60°C and a time range of 0-30 minutes, managed by a backlit central controller. Eurosonic® 3D operates with a Sweep Mode technology that distributes the ultrasonic waves evenly to achieve more consistent cleaning results and features a visual cycle count-down device which uses the LEDs that surround the front of the controller. The Eurosonic® 3D is fitted with two side handles to make it easier to move around.

- Large capacity stainless steel tank: 3 litres
- Overall dimensions: L 256 x H 220 x D 156 mm
- Weight: 3,3 kg
- Available accessories: perforated stainless steel basket, smoke-grey anti-deformation cover with non-drip system, 2 x 600 ml beaker holder and glass beakers.



Euroseal® 2001 Plus

Digital ultrasonic tank for dental clinics
Stainless steel ultrasonic tank with electromechanical set up. Equipped with Euronda Sweep Mode technology: working at different frequencies it guarantees that the ultrasonic waves are homogenously distributed all over the tank. The Energy is equipped with two side handles that facilitate the operations of emptying out of the ultrasonic tank. Control panel with easy to use on/off functions. Constant 60°C heating and working time adjustable from 0-30 minutes.

- Large capacity stainless steel tank: 3 litres
- Overall dimensions: L 256 x H 220 x D 156 mm
- Available accessories: plastic cover, perforated stainless steel basket, 2x600 ml beaker holder and 600 ml glass beaker.

Available accessories

- smoke-grey anti-deformation cover with non-drip system
- perforated stainless steel basket
- 2 x 600 ml beaker holder
- 600 ml glass beaker
- 1 x 600 ml beaker holder
- 3 x 600 ml beaker holder



Eurosonic® Micro

New compact ultrasonic tank for dentists

Electronic compact stainless steel ultrasonic tank with a single on-off function with a preset 10 minutes work cycle. Equipped with one piezoceramic transducer, Eurosonic® Micro guarantees a perfect cleaning, especially of little instruments (burs,etc.).

- Large capacity stainless steel tank: 0,5 litres
- Overall dimensions: L 235 x H 130 x D 170 mm
- Weight: 1 kg



X RAY ITEMS



Jolly box

Portable dark room for x-rays

Portable dark room for developing endoral x-rays.

Overall dimensions: L 320 x H 440 x D 290 mm



X-ray apron

Dental radiation front protection apron

Child X-ray apron with collar (Size: 480 x 530 mm)

- Lead protection 0,3 mm; weight 1,7 kg


Adult X-ray apron with collar (Size: 610 x 685 mm)

- Lead protection 0,3 mm, weight 2,8 kg

Operator X-ray apron (Size: 610 x 910 mm)

- Lead protection 0,3 mm; weight 3 kg

MONOART® DISPOSABLE LINE / PROTECTIVE GLASS

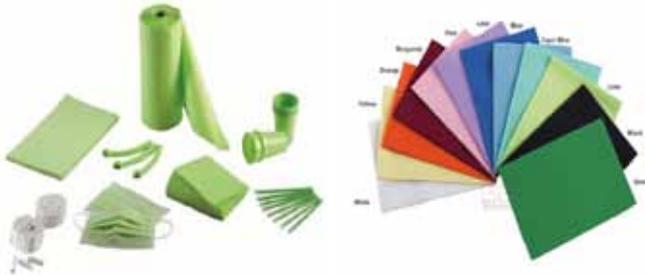

Monoart® products

Disposable products for 500 patients:

- 500 Monoart® Premium Towel (Dimensions: 330 x 450 mm) or 50 Monoart® Towel Up! (Dimensions: 330 x 450 mm)
- 500 Monoart® Plastic cup (Capacity: 166 cc / 200 cc)
- 500 Monoart® EM15 Saliva ejector (see below for more info.)
- 200 Monoart® Face mask Protection 3 (Pack: carton of 6 boxes)
- 250 Monoart® Tray paper (box of 250 pieces / 100 pieces)
- 160 Monoart® Apron PG30 (Size: roll of 80 aprons, size 610 x 530 mm; roll of 60 aprons, 810 x 530 mm)
- 10 Monoart® EM19 Aspirator tip (Tip: Ø 16 mm / length 128 mm)
- 600 g Monoart® Cotton roll (Size: Ø 8 mm L 37 mm, Ø 10 mm L 37 mm, Ø 12 mm L 37 mm)

Colours: light blue, green, yellow, orange, blue, pink, white, lilac, lime & etc.



Monoart® Infection Control

A safe barrier between you and infections

- 50 Monoart® Apron (Dimensions: 800 x 1100 mm, Colour: white)
- 50 Monoart® Face mask Protection 3 (Pack: carton of 6 boxes)
- Monoart® Bandana (Colours: 10 multicolours/ Box: single piece)
- Monoart® Ultra Light glasses (Weight: 29 g)
- Gown (Colour: white / Pack: carton of 25 pieces)
- Cut proof gloves (Colour: violet / Size: S, M / Pack: bag of 3 pairs)



Monoart® Cube Orange Glasses

PROTECTIVE GLASSES

Orange optically checked polycarbonate lens, providing a 180 degree view. Heavy duty lens, with antifog and antiscratch treatment, that protect against UV400 (sunlight) rays and UV25 (blue light) rays (polymerisation lamps). They can fit over any corrective glasses. Black and orange frame.

Weight: 37 g



Monoart® EM15 Saliva ejector

Saliva ejectors for dental clinic

Disposable flexible saliva ejector made of non-toxic transparent PVC. Integrated metallic wire for the excellent maintenance of the set angle. Rounded-shape fixed or removable tip, ideal for work on irritated mucous.

Size: diameter 6,5 mm

Medical Device: Class IIa

Length available: 125 or 150 mm